

NFT

Urine collection system

Needle-free, hygienic,
enclosed collection

Urine pre-
analytics
2.0



SARSTEDT

NFT

Needle Free Transfer

The first urine collection system from SARSTEDT for hygienic and needle-free sample collection

Hygienic. Practical. Safe.

How can collecting urine samples be made even safer, easier and more hygienic? The answer to this question can be found in the innovative products of the NFT system. As the first enclosed and needle-free urine collection system from SARSTEDT, the new product line combines the advantages of existing products in a single solution.


New standards for urine pre-analytics



- Available as an NFT urine cup, as well as an NFT 3-litre urine container for 24 hour urine collection
- Needle-free collection system
- Safe storage and transport
- Easy, hygienic urine transfer with the Urine-Monovette®
- No sharps container needed to dispose of the NFT products

Simply practical

The Urine-Monovette® enables hygienic urine collection and is also suitable as a primary tube for test strips and analysers and as a centrifuge tube for collecting sediment

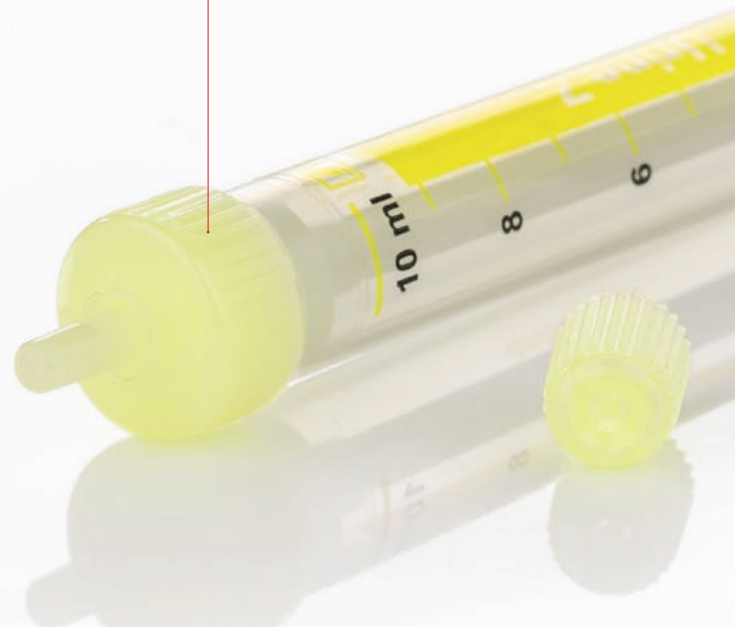


New

A needle is no longer needed with the new needle-free collection function. The Urine-Monovette® easily pierces the innovative NFT membrane.


Following an easy collection, the opening is immediately closed and hygienically sealed even after collecting several samples.

There is no risk of infection from a needlestick injury.

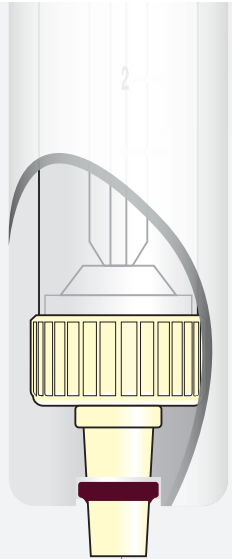

CONVINCING

- ✓ **Minimised risk of injury** and **increased user safety** for personnel and patient through needle-free sample collection
- ✓ **Savings on disposal**, as the components of the NFT system do not need to be disposed of separately in a Sharps container

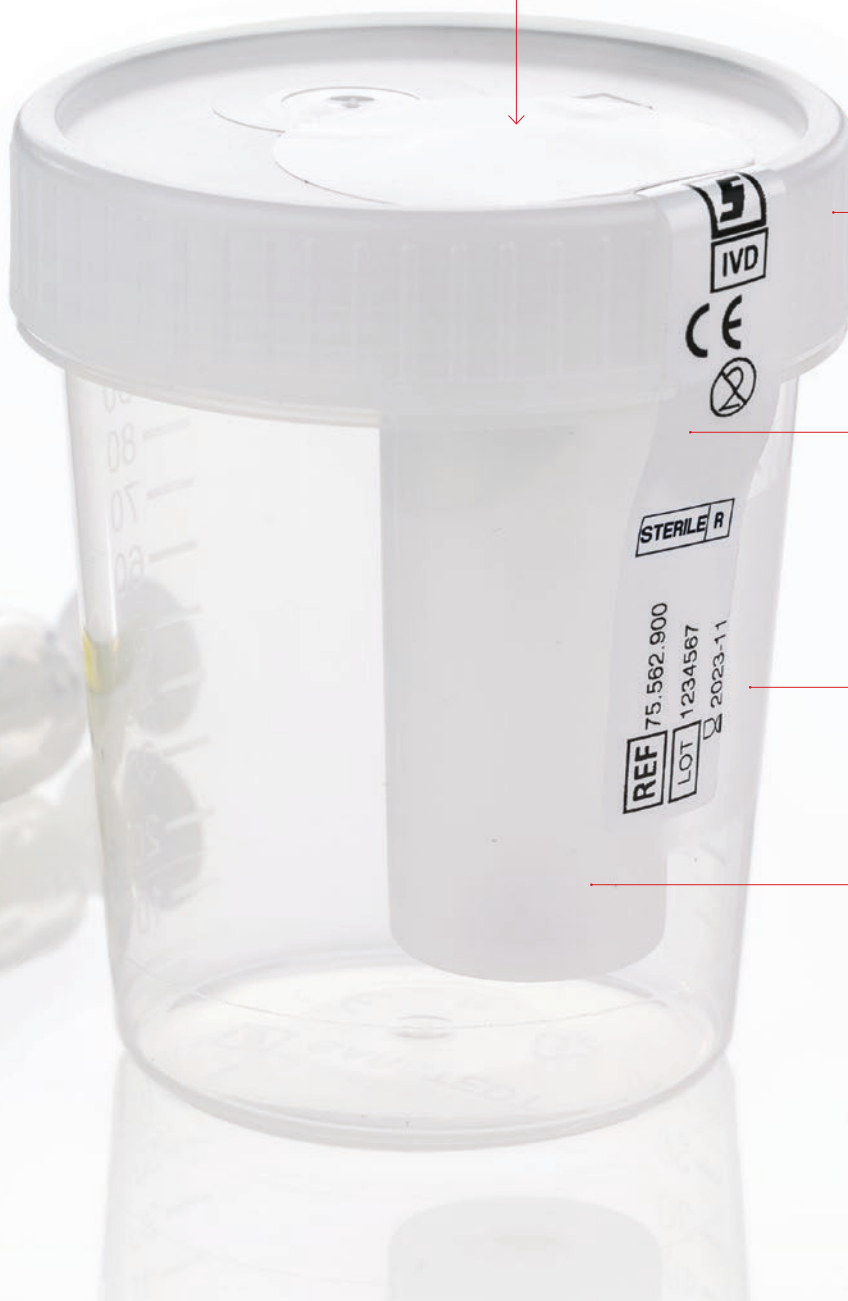

Liquid and odour tight
Screw cap with transfer unit

Sealed
If intact, the safety label indicates the internal sterility of the NFT urine cup

Clear
Highly transparent container of robust polypropylene (PP) for good visual assessment

Innovative
Needle-free transfer unit with pierceable NFT membrane

NFT

Urine collection system

The best from various systems

The NFT system with Urine-Monovette[®], NFT urine cup and NFT urine container ensures hygienic working conditions during the entire sample collection process:

- With the Urine-Monovette[®] and the enclosed suction tip, you can carry out the transfer process without a needle, even if the filling volume is very small. A separate transfer device as in vacuum systems is not needed.
- Extensive training is not required.
- With the Urine-Monovette[®], the aspiration principle allows the collection process to be repeated if necessary – there is no risk of loss of function through repeated use.
- The well-designed, needle-free collection function minimises the risk of injury for staff and patient.
- After collection, the products of the NFT system do not need to be disposed of separately in needle-proof containers – saving time and money compared with products with an integrated needle.

Well designed

The NFT urine container with a volume of 3,000 ml, ergonomic grip, large opening diameter and transparent inspection window with graduated scale in increments of 100 ml



Urine-Monovette® – hygienic work with just a few movements

- The Urine-Monovette® enables hygienic, needle-free sample collection from urine containers or urine drainage systems.
- A urine analysis with commercially available test strips can be carried out directly in the filled Urine-Monovette®.
- The Urine-Monovette® is suitable for transporting and shipping according to the P650 packaging instructions of the ADR and IATA. It can be used in centrifuges and analysers for collecting sediment.



User-friendly handling: NFT urine cup



Convenient for 24-hour urine collection: NFT urine container



The large-volume NFT urine container also ensures sample collection remains hygienic thanks to the enclosed collection function.



Watch the tutorial on the easy handling of the NFT products:
sarstedt.com/nft-tu-en

ORDER A FREE SAMPLE



sarstedt.com/nft-sa-en

SARSTEDT sets new standards with the NFT urine collection system

Only a clean and precise sample collection can ensure correct analysis results. The products of the NFT system from SARSTEDT set the standard for accurate pre-analytics suitable for routine use combined with the greatest safety for user and patient.

Urine pre-
analytics
2.0



NFT products do not need to be disposed of in Sharps containers.



SARSTEDT

NFT urine cup 100 ml

Order no.	Volume	L x Ø	Version	Purity	Packaging (inner/outer box)
75.562.900	100 ml	72 x 62 mm	With safety label and transfer unit	Sterile	200 / 200

NFT urine container 3.0 l

Order no.	Volume	L x W x H	Version	Purity	Packaging (inner/outer box)
77.575.900	3000 ml	125 x 130 x 246 mm	With enclosed transfer unit and patient label, opening diameter 80 mm	–	30 / 30

Urine-Monovette® Z - without additive

Order no.	Volume	L x Ø	Version (yellow)	Purity	Packaging (inner/outer box)
10.250.001	3,2 ml	75 x 13 mm	Paper label, writable	–	64 / 512
10.258	8,5 ml	92 x 15 mm	Plastic label, transparent	–	64 / 512
10.258.001	8,5 ml	92 x 15 mm	Paper label, writable	–	64 / 512
10.258.020	8,5 ml	92 x 15 mm	Plastic label, transparent, individually blister packed	Sterile	100 / 500
10.252	10 ml	102 x 15 mm	Plastic label, transparent	–	64 / 512
10.252.001	10 ml	102 x 15 mm	Paper label, writable	–	64 / 512
10.252.020	10 ml	102 x 15 mm	Plastic label, transparent, individually blister packed	Sterile	100 / 500
10.252.030	10 ml	102 x 15 mm	Plastic label, transparent, with light protection	–	64 / 512

Urine-Monovette® Boric acid - with stabiliser for microbiological tests

Order no.	Volume	L x Ø	Version (green)	Purity	Packaging (inner/outer box)
10.256*	3,2 ml	75 x 13 mm	Plastic label, transparent	–	64 / 512
10.256.020*	3,2 ml	75 x 13 mm	Plastic label, transparent, individually blister packed	Sterile	100 / 500
10.253	10 ml	102 x 15 mm	Plastic label, transparent	–	64 / 512
10.253.001	10 ml	102 x 15 mm	Paper label, writable	–	64 / 512
10.253.020	10 ml	102 x 15 mm	Plastic label, transparent, individually blister packed	Sterile	100 / 500
10.253.021	10 ml	102 x 15 mm	Paper label, writable, individually blister packed	Sterile	100 / 500

For more consumables for urine pre-analytics, see www.sarstedt.com

*Planned introduction in the third quarter of 2023

SARSTEDT AG & Co. KG

P.O. Box 12 20
D-51582 Nümbrecht

Phone: +49 2293 305 0
Fax: +49 2293 305 3992

Customer Service Germany
Phone +49 0800 0 83 305 0

export@sarstedt.com
www.sarstedt.com



Further information,
videos and sample order:
sarstedt.com/nft-en

If you have any questions:
We'd be happy to help!

Visit our website: www.sarstedt.com